



# **POSICIONADORES DE PIERNAS PARA MESAS QUIRÚRGICAS**

# INVENTARIO TÉCNICO

**Nombre:** POSICIONADORES DE PIERNAS PARA MESAS QUIRÚRGICAS


**Modelo:** S/M

**Marca:** STERIS

**Referencia:** 6509BF351

Cantidad	Rack	Columna	Tramo
7	L	1	1

Manual de instalación:

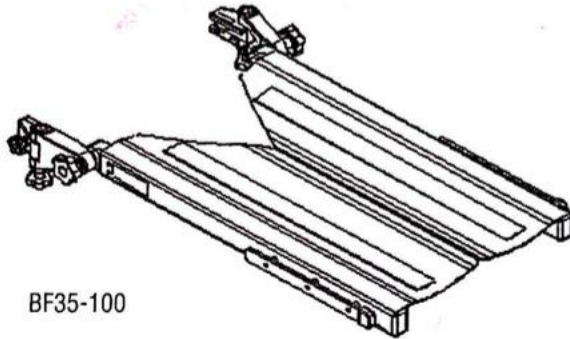

BF35-100

## **Amsco® Split Leg Positioner Operating Instructions**

### **Features**

- Lightweight
- Radiolucent Carbon Fiber Construction
- Compatible with Amsco® 2080, 3080, and 3085 Surgical Tables
- Contoured Leg Sections
- User Friendly
- Safe Patient Positioning

### **Procedures**

- Laparoscopies
- Low Anterior Resections
- Nissen Fundoplication
- Arthroscopies
- Lower Extremity Vascular and Orthopedic Procedures

**STERIS®**





Image 1



Image 2

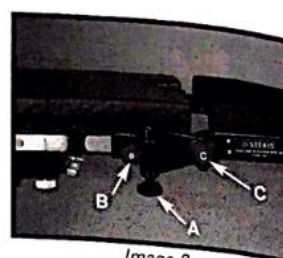


Image 3

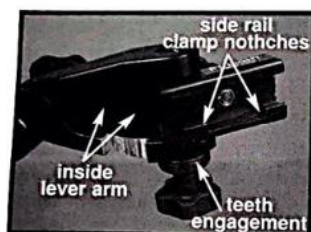


Image 4

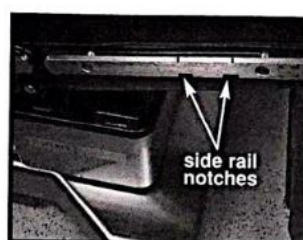


Image 5

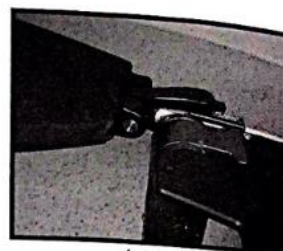


Image 6

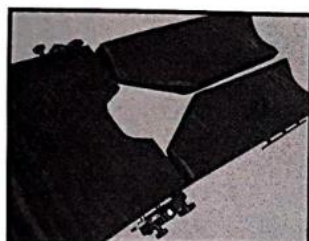


Image 7


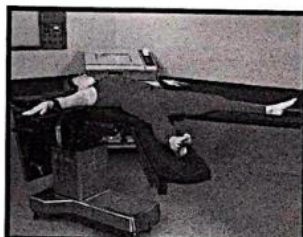


Image 8



Image 9

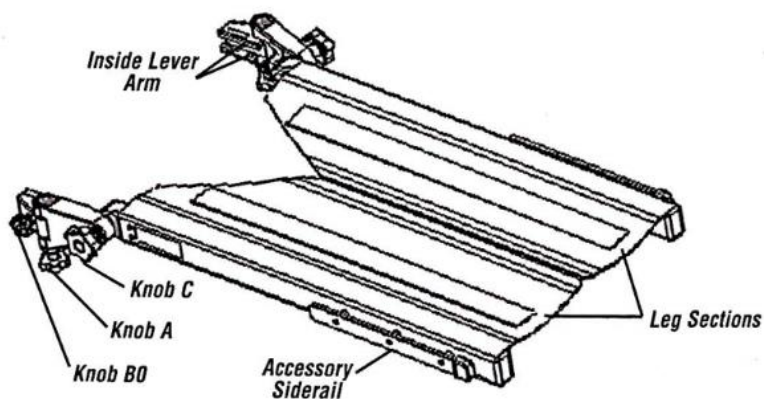

**Indications for use:** The Amsco Split Leg Positioner is a two section attachable accessory for Amsco surgical tables. It provides a safe system for positioning patients for laparoscopic procedures, as well as lower extremity orthopedic and vascular procedures. By utilizing the radiolucent Split Leg Positioner, both the surgeon and imaging equipment have excellent access to the operative site.

### **⚠ WARNING - SAFEGUARD PATIENT:**

Read and understand all instructions prior to using the Split Leg Positioner. If you need technical assistance or additional instructions, call STERIS Customer Service at 800-548-4873. Prior to using the Amsco Split Leg Positioner, assess the patient for any pre-existing conditions that might contraindicate the use of leg positioning devices.

## ATTACHING LEG SECTIONS

The Split Leg Positioner can be either attached onto the seat section of a 2080, 3080, or 3085 surgical table with the leg section lowered or, when the surgeon would prefer no obstruction from the surgical leg section, the accessory can be attached to the IA Extender on the 3080 and 3085 tables.

### STEP 1

**For attachment onto surgical table seat section**

Lower the surgical table leg section. *Refer to Image 1.*

**For attachment onto the 3080 or 3085 surgical table IA Extender**

Remove surgical table head section, attach IA Extender, ensure knobs are tightened, reverse the surgical table, and program the hand control accordingly. *Refer to Image 2.*

### STEP 2

Loosen the abduction knob (A) and set to position C (complete abduction 45°). *Refer to Images 3 and 4.*

### STEP 3

**For attachment onto surgical table seat section**

Align accessory side rail clamp notches with seat section side rail notches to attach leg section onto the surgical table. *Refer to Images 4 and 5.*

**For attachment onto the 3080 or 3085 surgical table IA Extender**

Slide the accessory side rail clamp onto the end of the IA Extender siderail to desired position.

**Note:** Located above the notches on the table siderail there are two lines which will help you align the clamp of the accessory with the notches in the siderail. *Refer to Image 5.*

### STEP 4

Once the accessory is attached to the side rail readjust abduction to position A (no abduction). *Refer to Image 6.*

### STEP 5

Slide accessory down the seat section or IA Extender side rail so that the end of the side rail is positioned inside the lever arm. *Refer to Images 4 and 6.*

### STEP 6

Once the Split Leg Positioner is in place tighten side rail knob (B) and ensure all other knobs are securely tightened and the leg section is parallel to the floor; ensure that teeth are engaged. *Refer to Images 4 and 7.*

### STEP 7

Repeat steps 2 through 6 for the opposite leg section.

## PATIENT POSITIONING

Once the patient is transferred onto the table and anesthetized, position patient's legs according to surgeon preference. Observe and protect any possible pressure or pinch points.

### TO RAISE OR LOWER LEG SECTION

While supporting the patient's leg loosen knob C, position leg to desired height, and tighten knob C, ensuring that teeth are engaged. *Refer to Images 3, 4, and 9.*


### TO ABDUCT LEG SECTION

Loosen knob A, abduct leg to desired position (A through C), and tighten knob A, ensuring that teeth are engaged. *Refer to Images 3, 4, and 8.*

### PATIENT LEG SUPPORT

When positioning patients for certain procedures each patient will require individual positioning according to weight and height. Once the legs are abducted to desired position it may be necessary to slide the accessory up the side rail towards the patient's head. This will allow more of the patient's leg to be supported by the accessory.



 **WARNING** Be aware of any pressure points and provide padding and support accordingly.

## COMPLETION OF PROCEDURES

### STEP 1

If the accessory was repositioned on the table side rail, loosen knob B and slide the Split Leg Positioner down to the original position with the side rail positioned just inside of the lever arm. *Refer to Images 3 and 6.*

### STEP 2

Both leg sections should be repositioned horizontally by loosening knob C and adjusting the height so that both legs are parallel to the floor. Ensure teeth are engaged and that knob C is retightened. *Refer to Images 4 and 8.*

### STEP 3

Loosen knob A to abduct the patient's legs back together to position A. *Refer to Image 7.*

## CLEANING INSTRUCTIONS

Clean each section, including the pads, with STERIS Coverage™ Spray HB, Coverage™ HB Concentrate, or Coverage Plus® prior to removing from the surgical table. Dry thoroughly when cleaning is complete.

## TO REMOVE LEG SECTIONS

Once the patient has been transferred to a stretcher, the leg sections can be removed.

## REMOVAL FROM SEAT SECTION SIDE RAIL

### STEP 1

Loosen knob A and abduct each section to the C position (complete abduction). *Refer to Image 4.*

### STEP 2

Loosen knob B and slide the leg section up the side rail, towards the head of the table, until the side rail clamp notches are aligned with the side rail notches. Lift the leg section off of side rail. *Refer to Images 4 and 5.*

## REMOVAL OF ATTACHMENT FROM 3080 AND 3085 IA EXTENSION

### STEP 1

Loosen knob B and slide accessory off the side rail.

## STORAGE

For storage, reposition each leg section to position A (no abduction) and tighten knobs A and C.

**STERIS®** STERIS Corporation  
5960 Heislley Road  
Mentor, OH 44060-1834 • USA  
440-354-2600 • 800-548-4873  
In Canada, 800-661-3937  
www.steris.com



**STERIS Has Offices Worldwide**  
Germany 49 7771 9332 0  
Hong Kong 852 2 563 3623  
Italy 39 0141 590429  
Japan 81 78 252 1901  
Korea 82 2 554 1661  
Latin America 1 305 442 8202  
Singapore 65 841 7677  
Spain 3491 658 5920  
United Kingdom 44 1 608 811822

## ITEMS

