



# **EXTENSIÓN DE MESA PARA MICROCIRUGÍA DE MANO Y BRAZO TRANSLÚCIDO A RX**



# INVENTARIO TÉCNICO

**Nombre:** EXTENSIÓN DE MESA PARA MICROCIROUGÍA DE MANO Y BRAZO TRANSLÚCIDO A RX

**Modelo:** S/M


**Marca:** STERIS

**Referencia:** 6509BF436

Cantidad	Rack	Columna	Tramo
10	I	1	2

## Manual de instalación:

### MICRO SURGICAL ARM AND HAND TABLES OPERATING INSTRUCTIONS

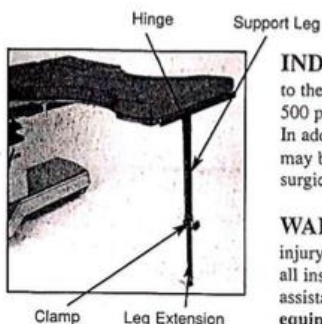

**INDICATIONS FOR USE** - The Micro Surgical Arm and Hand Table is an accessory that attaches to the surgical table side rail with accessory clamps. It is intended to support the operative arm of up to 500 pound patients, provide a stable and radiolucent surgical platform, and provide elbow rest for user(s). In addition, the side rail is designed to support traction accessories. (Note: Surgical table load capacities may be less. **Never** overload a surgical table.) The device is intended for mounting on the side rail of a surgical table only.

**WARNING - SAFEGUARD PATIENT:** Cracked or splintered surfaces may cause injury. Inspect table top prior to use. Damaged table top must be replaced. Read and understand all instructions prior to using the Micro Surgical Arm and Hand Table. If you need technical assistance or additional instructions, call STERIS Customer Service at 800-548-4873. **Do not use equipment if worn, damaged, or pieces are missing.**

#### Attachment/Set-up Instructions

- Step 1. - Attach two STERIS Clark Sockets or STERIS Accessory Clamps onto the operative side of the surgical table side rails. Place these approximately 13 inches apart from each other and in the general location of the patient's shoulder.
- Step 2. - Insert the posts of the Arm and Hand Table into each clamp.
- Step 3. - Slide the table and clamps to the desired location on the side rails then lock the arm and hand table into the clamp and lock onto the side rail.

**WARNING** - To prevent operator injury, keep hands and fingers away from the hinge of the support leg.

- Step 4. - Fully extend the support leg out until it reaches its operating position. Lower the leg extension to the floor and lock the leg extension to effectively support the weight.
- Step 5. - To adjust the height of the Arm and Hand Table, loosen the accessory clamps and raise or lower the table to the desired location, then lock the clamps. Readjust the support leg height as needed to effectively support the weight.

**WARNING** - To prevent patient or operator injury, ensure the clamps are locked into place and the Arm and Hand Table posts are safely secured into the clamp.

#### Removal Instructions

**WARNING** - To prevent operator injury, keep hands and fingers away from the hinge of the support leg.

- Step 1. - Loosen the accessory clamps and remove the Arm and Hand Table by lifting up on the board.
- Step 2. - Retract and lock the leg extension then press the support leg against the lock underneath the Arm and Hand Table for storage.
- Step 3. - Remove the accessory clamps from the side rails.

#### Cleaning Instructions

Clean the Arm and Hand Table and pad surface with STERIS Coverage™ Spray HB, Coverage™ HB Concentrate, or Coverage Plus®. Dry assembly thoroughly when cleaning is complete.

Publication ID #41913-000110 ©2008 by STERIS Corporation. All rights reserved.	
STERIS Offices Worldwide	
Canada	800 561 3037
Germany	49 2333 5999 400
Hong Kong	852 2 563 3623
Italy	39 0141 58 54 29
Japan	81 19 262 1561
Korea	82 2 554 1661
Latin America	1 302 442 8202
Singapore	65 641 1077
Spain	34 918 563 920
United Kingdom	44 1 858 911822
Infection and Contamination Prevention, Worldwide	



STERIS Corporation  
5960 Halsey Road  
Mentor, OH 44060-1834 • USA  
440-354-2600 • 800-548-4873  
www.steris.com