



ACCESORIOS PARA MESA CLÍNICA

INVENTARIO TÉCNICO

Nombre: ACCESORIOS PARA MESA CLÍNICA

Modelo: EMAX 300 Y ASC 2000

Marca: STERIS

Referencia: BF-643

Cantidad	Rack	Columna	Tramo
2	K	2	1

Manual de instalación:

STERIS®



OPERATOR MANUAL

Bariatric Table Extensions

(03/31/06)

P150830-480

A WORD FROM STERIS CORPORATION

This manual contains important information on proper use and maintenance of the Bariatric Table Extensions. **All personnel involved in the use and maintenance of this equipment must carefully review and comply with the warnings, cautions, and instructions contained in this manual.** These instructions are important to protect the health and safety of patients and personnel operating these extensions. This manual should be retained in a conveniently accessible area for quick reference.

STERIS carries a complete line of surgical table accessories. STERIS representatives will gladly review these with you.

Advisory

A listing of the *Safety Precautions* to be observed when operating and servicing this equipment can be found in *Section 1* of this manual. Do not operate the equipment until you have become familiar with this information.

Any alteration of this equipment not authorized or performed by STERIS that could affect unit operation will void the warranty and could violate federal, state, and local regulations.

Indications for Use

The Bariatric Table Extensions are an accessory, attachable to the head, back, seat, and/or legs sections of Amsco® 3085 SP, 3080, and Cmax™ Surgical Tables **only**. Users must verify the patient weight capacity of the surgical table with this accessory applied.

The extensions (when properly applied to each side of the operating table) increase the table width by a total of eight inches (203 mm) and provide a stable platform for patient support and the attachment of many Amsco 3080/3085 SP and Cmax Surgical Table accessories.

These extensions are available in various sets as follows:

- Full Set with either 2" or 4" pads
- Seat/Back Set with either 2" or 4" pads
- Leg Only with either 2" or 4" pads
- Head Only with either 2" or 4" pads

Each extension is clearly marked to specify the type of extension (head, back, seat, leg) and the table designed for (Amsco 3080/3085 SP or Cmax). Also, the extension pads for the Cmax Surgical Table are identified with blue attachment strips.

TABLE OF CONTENTS






Section	Title	Page
	A WORD FROM STERIS CORPORATION	i
1	SAFETY PRECAUTIONS	1-1
2	OPERATING INSTRUCTIONS	2-1
	2.1 Before Using Bariatric Table Extensions	2-1
	2.2 Attaching Bariatric Table Extensions	2-3
	2.3 Transferring Patient to Table	2-4
	2.4 Removing Bariatric Table Extensions	2-4
3	ROUTINE MAINTENANCE	3-1
	3.1 Cleaning Bariatric Table Extensions	3-1
	3.2 Replacement Parts	3-1
4	DISPOSAL HAZARDS	4-1

SAFETY PRECAUTIONS



1

The following list of *Safety Precautions* **must** be observed when operating and servicing the Amsco® 3080/3085 SP or Cmax™ Surgical Table Bariatric Table Extensions. WARNING indicates the potential for personal injury and CAUTION indicates the potential for damage to equipment. For emphasis, certain *Safety Precautions* are repeated throughout the manual. It is important to review ALL *Safety Precautions* before operating this equipment.


WARNING – PERSONAL INJURY AND/OR EQUIPMENT DAMAGE HAZARD:

-  To prevent patient and/or user injury and/or equipment damage, read and understand all instructions presented in this manual before using the Bariatric Table Extensions. Follow each step in the order presented in these instructions. If you need technical assistance or additional instructions, call STERIS.
-  To prevent patient and/or user injury and/or equipment damage, check the Bariatric Table Extensions and surgical table side rails for possible damage or wear prior to each use. Do not use the Bariatric Table Extensions if any damage is apparent, if parts are missing, or if it does not operate as expected. Cracked or splintered surfaces may cause injury.
-  To prevent possible patient or user injury or instability of the table, apply extensions to BOTH sides of the table section. DO NOT use one extension per section.
-  To prevent patient and/or user injury and/or equipment damage, apply Bariatric Table Extensions **only** to Amsco 3080/3085 SP or Cmax surgical table side rails. Amsco 3080/3085 SP and Cmax Bariatric Table Extensions ARE NOT interchangeable.
-  To prevent patient and/or user injury and/or equipment damage, assure Bariatric Table Extensions attaching clamps fully contact table side rails and are securely tightened in place. Test locking to ensure no extension movement when lifted or pushed.



WARNING – TIPPING HAZARD:

-  Do not place a patient on the surgical table unless the table floor locks are engaged.
-  Do not release the surgical table floor locks while the patient is on the table.

WARNING – INFECTION HAZARD:

-  When cleaning the Bariatric Table Extensions, follow the cleaning product manufacturer's instructions and precautions for the cleaning agent being used.

CAUTION – POSSIBLE EQUIPMENT DAMAGE:

-  Handle the Bariatric Table Extensions with care. Avoid dropping the equipment to prevent damage. Do not use if damaged.
-  To prevent patient and/or user injury and/or equipment damage, DO NOT attach a Bariatric Table Extension to another Bariatric Table Extension.

OPERATING INSTRUCTIONS

2

2.1 Before Using Bariatric Table Extensions

⚠ WARNING – PERSONAL INJURY AND/OR EQUIPMENT DAMAGE:

- To prevent patient and/or user injury and/or equipment damage, read and understand all instructions presented in this manual before using the Bariatric Table Extensions. Follow each step in the order presented in these instructions. If you need technical assistance or additional instructions, call STERIS.
- To prevent patient and/or user injury and/or equipment damage, apply Bariatric Table Extensions only to Amsco 3080/3085 SP or Cmax surgical table side rails. Amsco 3080/3085 SP and Cmax Bariatric Table Extensions ARE NOT interchangeable.

The following instructions are intended as a guide for using the Bariatric Table Extensions.

NOTE: Refer to the operator manual for your Amsco® 3085 SP or 3080 Series or Cmax™ Surgical Table for specific instructions on operating the surgical table.

The Bariatric Table Extensions (see Fig. 2-1) are designed as a user installed modular system on Amsco 3080/3085 SP and Cmax Surgical Tables. When properly applied to each side of the operating table, the extensions increase the original table width by a total of eight inches (203 mm). Bariatric Table Extensions provide a stable platform for patient support and the attachment of many 3080/3085 SP and Cmax Surgical Table accessories. The Bariatric Table Extensions system is available in various sets for facility flexibility as follows:

1. The complete (full set) Bariatric Table Extensions system is comprised of eight pieces identified as: *Head Section, Back Section, Seat Section, and Leg Section*. The extensions fit either side of the identified table section.

NOTE: The Bariatric Extensions are unique to each table section and are not interchangeable with other table sections. The Amsco 3080/3085 SP and Cmax Bariatric Table Extensions ARE NOT interchangeable.

2. The Seat/Back Set is comprised of four pieces identified as: *Back Section and Seat Section*. The extensions fit either side of the identified table section.
3. The Leg Only Set is comprised of two pieces identified as *Leg Section*. The extensions fit either side of the identified table section.

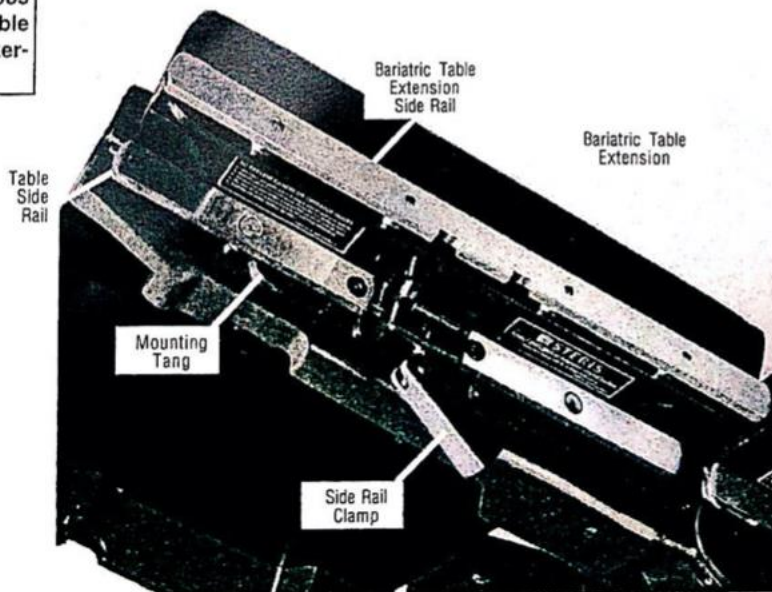


Figure 2-1. Amsco 3080/3085 SP Bariatric Table Extension Components

2.2 Attaching Bariatric Table Extensions

WARNING – PERSONAL INJURY AND/OR EQUIPMENT DAMAGE:

- To prevent patient and/or user injury and/or equipment damage, check the Bariatric Table Extensions and surgical table side rails for possible damage or wear prior to each use. Do not use the Bariatric Table Extensions if any damage is apparent, if parts are missing, or if it does not operate as expected. Cracked or splintered surfaces may cause injury.
- To prevent patient and/or user injury and/or equipment damage, apply Bariatric Table Extensions only to Amsco 3080/3085 SP or Cmax surgical table side rails. Amsco 3080/3085 SP and Cmax Bariatric Table Extensions ARE NOT interchangeable.
- To prevent possible patient or user injury or instability of the table, apply extensions to BOTH sides of the table section. DO NOT use one extension per section.
- To prevent patient and/or user injury and/or equipment damage, assure Bariatric Table Extensions attaching clamps fully contact table side rails and are securely tightened in place. Test locking to ensure no extension movement when lifted or pushed.

CAUTION – POSSIBLE EQUIPMENT HAZARD:

- Handle the Bariatric Table Extensions with care. Avoid dropping the equipment to prevent damage. Do not use if damaged.
- To prevent patient and/or user injury and/or equipment damage, DO NOT attach a Bariatric Table Extension to another Bariatric Table Extension.

Complete each step in the order presented below and attach the Bariatric Table Extensions to an Amsco 3080/3085 SP or Cmax Surgical Table as follows:

1. Check bottom of Bariatric Table Extensions for identifier labels (i.e., Head, Back, Seat, Leg Section) to ensure proper attachment to correct table section. **Verify each surgical table side rail contains no other accessories before attaching extension.**

2. Ensure extension side rail clamp (see Fig. 2-1) jaw is fully open by rotating clamp locking handle in a counter-clockwise motion (as viewed from bottom of extension).

3. Position correct extension above table side rail and insert extension mounting tang between table top and side rail (see Fig. 2-2). Lower extension into place, resting on top of table side rail.

NOTE: The extension should fit evenly over the surgical table side rail and appear to be a continuation of the tabletop itself.

4. Tighten extension side rail clamp by rotating clamp locking handle clockwise (as viewed from below).

5. Install appropriate pads to top of Bariatric Table Extension.

NOTE: For other table accessories, note the following:

- 1) Split Leg Positioners can be attached to the Bariatric Table Extensions for mid-line access to the patient's abdomen.
- 2) The new Bariatric Foot Extension is designed for use with the Bariatric Table Extensions (replacing the narrower Foot Extensions).
- 3) Most accessories can still be used on the Bariatric Table Extensions except for accessories that must be attached simultaneously to both of the extended side rails.

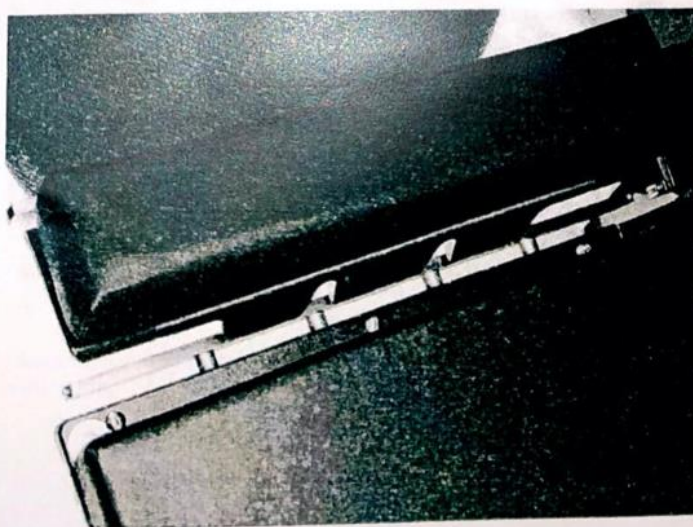


Figure 2-2. Typical Bariatric Table Extension Attachment – Amsco 3080 Series Surgical Table Shown

ROUTINE MAINTENANCE

3

3.1 Cleaning Bariatric Table Extensions



WARNING—INFECTION HAZARD: When cleaning the Bariatric Table Extensions, follow the cleaning product manufacturer's instructions and precautions for the cleaning agent being used.

1. Clean the Bariatric Table Extensions after each use.
2. STERIS recommends Coverage Plus NPD® One-Step Cleaner Disinfectant and either Coverage® Spray HB Plus Ready-to-Use One-Step Disinfectant or Coverage® Spray TB Plus Ready-to-Use Disinfectant Cleaner.* See Table 3-1 for a list of recommended cleaning products.
3. When cleaning, thoroughly read the cleaning product's instructions for use and follow all directions and precautions on the label. Pay particular attention to the following areas where fluid migration could occur:
 - handles
 - side rail clamps
4. After cleaning, ensure Bariatric Table Extensions and pads are completely dry before reinstalling, to prevent damage.

Table 3-1. Recommended Cleaning Products*

Product Name	Product Use
Coverage® Spray TB Plus Ready-to-Use Disinfectant Cleaner	Hospital-grade quaternary-based disinfectant spray
Coverage® Spray HB Plus Ready-to-Use Disinfectant	EPA-registered Hepatitis-B-effective quat
T.B.Q.® Hard Surface Disinfectant	Detergent-based germicide
Coverage Plus NPD® One-Step Cleaner Disinfectant	One-step disinfectant cleaner
Disinfectant Cleaner	General cleaner/disinfectant formulated with quaternary ammonium compounds and nonionic detergents
Germicidal Surface Wipes	Surface disinfectant

*Available from STERIS.

3.2 Replacement Parts

Please contact STERIS customer service (1-800-333-8828) for listing of replacement parts.

DISPOSAL HAZARDS

4

As of the date of this manual, materials (e.g., Herculite® Fusion®¹, polyurethane foam, alloy steel, stainless steel, aluminum, and adhesives) used in construction of the Bariatric Table Extensions and pads are not known to require special disposal methods. When disposal of the Extensions and/or pads become necessary, STERIS recommends the customer review the federal, state, and local regulations in your area to ensure proper disposal.

¹Herculite and Fusion are registered trademarks of Herculite, Inc.

**Protect your STERIS equipment with
cost-effective extended services agreements**

The best way to prevent costly downtime due to equipment malfunction is with regularly scheduled maintenance performed by qualified technicians trained in the latest technology. STERIS offers annual maintenance agreements to give your capital equipment planned maintenance that will help correct little problems before they become big ones. STERIS Engineering Service combines the precise maintenance program and factory-trained technicians to assure you of maximum productivity.

Our STERIS service technicians thoroughly inspect, clean, adjust and provide all necessary maintenance to keep your equipment performing according to factory specifications, all at an established economical rate that you can plan for.

We have more than 5,000 customers who are benefiting from STERIS maintenance agreements. Why not join them? Obtain complete details by calling 1-800-333-8828, or writing to:

STERIS®



STERIS Corporation
5960 Heisley Road
Mentor, OH 44060-1834 • USA
440-354-2600 • 800-JIT-4-USE (800-548-4873)